

laboratory, either manually by puncturing the septum with a small-gauge needle or automatically by the sample introduction system, just prior to analysis.

6.1.2 High concentration soil samples collected without a preservative

When high concentration samples are collected without a preservative, a variety of sample containers may be employed, including 60-mL glass vials with septum seals (see Sec. 4.4).

6.1.3 High concentration soil samples collected and preserved in the field

The following steps apply to the preparation of vials used in the collection of high concentration soil samples to be preserved in the field with methanol and analyzed by the aqueous purge-and-trap equipment described in Method 5030.

6.1.3.1 Add 10 mL of methanol to each vial.

6.1.3.2 Seal the vial with the screw-cap and septum seal.

6.1.3.3 Affix a label to each vial. This eliminates the need to label the vials in the field and assures that the tare weight of the vial includes the label. (The weight of any markings added to the label in the field is negligible).

6.1.3.4 Weigh the prepared vial to the nearest 0.01 g, record the tare weight, and write it on the label.

NOTE: Vials containing methanol should be weighed a second time on the day that they are to be used. Vials found to have lost methanol (reduction in weight of >0.01 g) should not be used for sample collection.

6.1.3.5 Surrogates, internal standards and matrix spikes (if applicable) should be added to the sample after it is returned to the laboratory and prior to analysis.

6.1.4 Oily waste samples

When oily waste samples are known to be soluble in methanol or PEG, sample vials may be prepared as described in Sec. 6.1.3, using the appropriate solvent. However, when the solubility of the waste is unknown, the sample should be collected without the use of a preservative, in a vial such as that described in Sec. 6.1.2.

6.2 Sample collection

Collect the sample according to the procedures outlined in the sampling plan. As with any sampling procedure for volatiles, care must be taken to minimize the disturbance of the sample in order to minimize the loss of the volatile components. Several techniques may be used to transfer a sample to the relatively narrow opening of the low concentration soil vial. These include devices such as the EnCore™ sampler, the Purge-and-Trap Soil Sampler™, and a cut plastic syringe. Always wear gloves whenever handling the tared sample vials.

6.2.1 Low concentration soil samples

6.2.1.1 Using an appropriate sample collection device, collect approximately 5 g of sample as soon as possible after the surface of the soil or other solid material has been exposed to the atmosphere: generally within a few minutes at most. Carefully wipe the exterior of the sample collection device with a clean cloth or towel.

6.2.1.2 Using the sample collection device, add about 5 g (2 - 3 cm) of soil to the sample vial containing the preservative solution. Quickly brush any soil off the vial threads and immediately seal the vial with the septum and screw-cap. Store samples on ice at 4°C.

NOTE: Soil samples that contain carbonate minerals (either from natural sources or applied as an amendment) may effervesce upon contact with the acidic preservative solution in the low concentration sample vial. If the amount of gas generated is very small (i.e., several mL), any loss of volatiles as a result of such effervescence may be minimal if the vial is sealed quickly. However, if larger amounts of gas are generated, not only may the sample lose a significant amount of analyte, but the gas pressure may shatter the vial if the sample vial is sealed. Therefore, when samples are known or suspected to contain high levels of carbonates, a test sample should be collected, added to a vial, and checked for effervescence. If a rapid or vigorous reaction occurs, discard the sample and collect low concentration samples in vials that do not contain the preservative solution.

6.2.1.3 When practical, use a portable balance to weigh the sealed vial containing the sample to ensure that 5.0 ± 0.5 g of sample were added. The balance should be calibrated in the field using an appropriate weight for the sample containers employed (Sec. 4.5.5). Record the weight of the sealed vial containing the sample to the nearest 0.01 g.

6.2.1.4 Alternatively, collect several trial samples with plastic syringes. Weigh each trial sample and note the length of the soil column in the syringe. Use these data to determine the length of soil in the syringe that corresponds to 5.0 ± 0.5 g. Discard each trial sample.

6.2.1.5 As with the collection of aqueous samples for volatiles, collect at least two replicate samples. This will allow the laboratory an additional sample for reanalysis. The second sample should be taken from the same soil stratum or the same section of the solid waste being sampled, and within close proximity to the location from which the original sample was collected.

6.2.1.6 In addition, since the soil vial cannot be opened without compromising the integrity of the sample, at least one additional aliquot of sample must be collected for screening, dry weight determination, and high concentration analysis (if necessary). This third aliquot may be collected in a 60-mL glass vial or a third 40-mL soil sample vial. However, this third vial must *not* contain the sample preservative solution, as an aliquot will be used to determine dry weight. If high concentration samples are collected in vials containing methanol, then two additional aliquots should be collected, one for high concentration analysis collected in a vial containing methanol, and another for the dry weight determination in a vial without either methanol or the low concentration aqueous preservative solution.

6.2.1.7 If samples are known or expected to contain target analytes over a wide range of concentrations, thereby requiring the analyses of multiple sample aliquots, it may be advisable and practical to take an additional sample aliquot in a low concentration soil vial containing the preservative, but collecting only 1-2 g instead of the 5 g collected in Sec. 6.2.1.1. This aliquot may be used for those analytes that exceed the instrument calibration range in the 5-g analysis.

6.2.1.8 The EnCore™ sampler has not been thoroughly evaluated by EPA as a sample storage device. While preliminary results indicate that storage in the EnCore™ device may be appropriate for up to 48 hours, samples collected in this device should be transferred to the soil sample vials as soon as possible, or analyzed within 48 hours.

6.2.1.9 The collection of low concentration soil samples in vials that contain methanol is not appropriate for samples analyzed with the closed-system purge-and-trap equipment described in this method (see Sec. 6.2.2).

6.2.2 High concentration soil samples preserved in the field

The collection of soil samples in vials that contain methanol has been suggested by some as a combined preservation and extraction procedure. However, this procedure is not appropriate for use with the low concentration soil procedure described in this method.

NOTE: The use of methanol preservation has not been formally evaluated by EPA and analysts must be aware of two potential problems. First, the use of methanol as a preservative and extraction solvent introduces a significant dilution factor that will raise the method quantitation limit beyond the operating range of the low concentration direct purge-and-trap procedure (0.5-200 µg/kg). The exact dilution factor will depend on the masses of solvent and sample, but generally exceeds 1000, and may make it difficult to demonstrate compliance with regulatory limits or action levels for some analytes. Because the analytes of interest are volatile, the methanol extract cannot be concentrated to overcome the dilution problem. Thus, for samples of unknown composition, it may still be necessary to collect an aliquot for analysis by this closed-system procedure and another aliquot preserved in methanol and analyzed by other procedures. The second problem is that the addition of methanol to the sample is likely to cause the sample to fail the ignitability characteristic, thereby making the unused sample volume a hazardous waste.

6.2.2.1 When samples are known to contain volatiles at concentrations high enough that the dilution factor will not preclude obtaining results within the calibration range of the appropriate determinative method, a sample may be collected and immediately placed in a sample vial containing purge-and-trap grade methanol.

6.2.2.2 Using an appropriate sample collection device, collect approximately 5 g of sample as soon as possible after the surface of the soil or other solid material has been exposed to the atmosphere: generally within a few minutes at most. Carefully wipe the exterior of the sample collection device with a clean cloth or towel.

6.2.2.3 Using the sample collection device, add about 5 g (2 - 3 cm) of soil to the vial containing 10 mL of methanol. Quickly brush any soil off the vial threads and immediately seal the vial with the septum and screw-cap. Store samples on ice at 4°C.

6.2.2.4 When practical, use a portable balance to weigh the sealed vial containing the sample to ensure that 5.0 ± 0.5 g of sample were added. The balance should be calibrated in the field using an appropriate weight for the sample containers employed (Sec. 4.5.5). Record the weight of the sealed vial containing the sample to the nearest 0.01 g.

6.2.2.5 Alternatively, collect several trial samples with plastic syringes. Weigh each trial sample and note the length of the soil column in the syringe. Use these data to determine the length of soil in the syringe that corresponds to 5.0 ± 0.5 g. Discard each trial sample.

6.2.2.6 Other sample weights and volumes of methanol may be employed, provided that the analyst can demonstrate that the sensitivity of the overall analytical procedure is appropriate for the intended application.

6.2.2.7 The collection of at least one additional sample aliquot is required for the determination of the dry weight, as described in Sec. 6.2.1.6. Samples collected in methanol should be shipped as described in Sec. 6.3, and must be clearly labeled as containing methanol, so that the samples are not analyzed using the closed-system purge-and-trap equipment described in this procedure.

6.2.3 High concentration soil sample not preserved in the field

The collection of high concentration soil samples that are not preserved in the field generally follows similar procedures as for the other types of samples described in Secs. 6.2.1 and 6.2.2, with the obvious exception that the sample vials contain neither the aqueous preservative solution nor methanol. However, when field preservation is not employed, it is better to collect a larger volume sample, filling the sample container as full as practical in order to minimize the headspace. Such collection procedures generally do not require the collection of a separate aliquot for dry weight determination, but it may be advisable to collect a second sample aliquot for screening purposes, in order to minimize the loss of volatiles in either aliquot.

6.2.4 Oily waste samples

The collection procedures for oily samples depend on knowledge of the waste and its solubility in methanol or other solvents.

6.2.4.1 When an oily waste is known to be soluble in methanol or PEG, the sample may be collected in a vial containing such a solvent (see Sec. 6.1.4), using procedures similar to those described in Sec. 6.2.2.

6.2.4.2 When the solubility of the oily waste is not known, the sample should either be collected in a vial without a preservative, as described in Sec. 6.2.3, or the solubility of a trial sample should be tested in the field, using a vial containing solvent. If the trial sample is soluble in the solvent, then collect the oily waste sample as described in Sec. 6.2.2. Otherwise, collect an unpreserved sample as described in Sec. 6.2.3.

6.3 Sample handling and shipment

All samples for volatiles analysis should be cooled to approximately 4°C, packed in appropriate containers, and shipped to the laboratory on ice, as described in the sampling plan.

6.4 Sample storage

6.4.1 Once in the laboratory, store samples at 4°C until analysis. The sample storage area should be free of organic solvent vapors.

6.4.2 All samples should be analyzed as soon as practical, and within the designated holding time from collection. Samples not analyzed within the designated holding time must be noted and the data are considered minimum values.

6.4.3 When the low concentration samples are strongly alkaline or highly calcareous in nature, the sodium bisulfate preservative solution may not be strong enough to reduce the pH of the soil/water solution to below 2. Therefore, when low concentration soils to be sampled are known or suspected to be strongly alkaline or highly calcareous, additional steps may be required to preserve the samples. Such steps include: addition of larger amounts of the sodium bisulfate preservative to non-calcareous samples, storage of low concentration samples at -10°C (taking care not to fill the vials so full that the expansion of the water in the vial breaks the vial), or significantly reducing the maximum holding time for low concentration soil samples. Whichever steps are employed, they should be clearly described in the sampling and QA project plans and distributed to both the field and laboratory personnel. See Sec. 6.2.1.2 for additional information.

7.0 PROCEDURE

This section describes procedures for sample screening, the low concentration soil method, the high concentration soil method, and the procedure for oily waste samples. High concentration samples are to be introduced into the GC system using Method 5030. Oily waste samples are to be introduced into the GC system using Method 5030 if they are soluble in a water-miscible solvent, or using Method 3585 if they are not.

7.1 Sample screening

7.1.1 It is highly recommended that all samples be screened prior to the purge-and-trap GC or GC/MS analysis. Samples may contain higher than expected quantities of purgeable organics that will contaminate the purge-and-trap system, thereby requiring extensive cleanup and instrument maintenance. The screening data are used to determine which is the appropriate sample preparation procedure for the particular sample, the low concentration closed-system direct purge-and-trap method (Sec. 7.2), the high concentration (methanol extraction) method (Sec. 7.3), or the nonaqueous liquid (oily waste) methanol or PEG dilution procedure (Sec. 7.4).

7.1.2 The analyst may employ any appropriate screening technique. Two suggested screening techniques employing SW-846 methods are:

7.1.2.1 Automated headspace (Method 5021) using a gas chromatograph (GC) equipped with a photoionization detector (PID) and an electrolytic conductivity detector (HECD) in series, or,

7.1.2.2 Extraction of the sample with hexadecane (Method 3820) and analysis of the extract on a GC equipped with a FID and/or an ECD.

7.1.3 The analyst may inject a calibration standard containing the analytes of interest at a concentration equivalent to the upper limit of the calibration range of the low concentration soil method. The results from this standard may be used to determine when the screening results approach the upper limit of the low concentration soil method. There are no linearity or other performance criteria associated with the injection of such a standard, and other approaches may be employed to estimate sample concentrations.

7.1.4 Use the low concentration closed-system purge-and-trap method (Sec. 7.2) if the estimated concentration from the screening procedure falls within the calibration range of the selected determinative method. If the concentration exceeds the calibration range of the low concentration soil method, then use either the high concentration soil method (Sec. 7.3), or the oily waste method (Sec. 7.4).

7.2 Low concentration soil method (Approximate concentration range of 0.5 to 200 µg/kg - the concentration range is dependent upon the determinative method and the sensitivity of each analyte.)

7.2.1 Initial calibration

Prior to using this introduction technique for any GC or GC/MS method, the system must be calibrated. General calibration procedures are discussed in Method 8000, while the determinative methods and Method 5000 provide specific information on calibration and preparation of standards. Normally, external standard calibration is preferred for the GC methods (non-MS detection) because of possible interference problems with internal standards. If interferences are not a problem, or when a GC/MS method is used, internal standard calibration may be employed.

7.2.1.1 Assemble a purge-and-trap device that meets the specification in Sec. 4.2 and that is connected to a gas chromatograph or a gas chromatograph/mass spectrometer system.

7.2.1.2 Before initial use, a Carbowack/Carbosieve trap should be conditioned overnight at 245°C by backflushing with an inert gas flow of at least 20 mL/minute. If other trapping materials are substituted for the Carbowack/Carbosieve, follow the manufacturers recommendations for conditioning. Vent the trap effluent to the hood, not to the analytical column. Prior to daily use, the trap should be conditioned for 10 minutes at 245°C with backflushing. The trap may be vented to the analytical column during daily conditioning; however, the column must be run through the temperature program prior to analysis of samples.

7.2.1.3 If the standard trap in Sec. 4.2.2.2 is employed, prior to initial use, the trap should be conditioned overnight at 180°C by backflushing with an inert gas flow of at least 20 mL/min, or according to the manufacturer's recommendations. Vent the trap effluent to the hood, not to the analytical column. Prior to daily use, the trap should be conditioned for 10 min at 180°C with backflushing. The trap may be vented to the analytical column during daily conditioning; however, the column must be run through the temperature program prior to analysis of samples.

7.2.1.4 Establish the purge-and-trap instrument operating conditions. Adjust the instrument to inject 5 mL of water, to heat the sample to 40°C, and to hold the sample at 40°C for 1.5 minutes before commencing the purge process, or as recommended by the instrument manufacturer.

7.2.1.5 Prepare a minimum of five initial calibration standards containing all the analytes of interest and surrogates, as described in Method 8000, and following the instrument manufacturer's instructions. The calibration standards are prepared in organic-free reagent water. The volume of organic-free reagent water used for calibration must be the same volume used for sample analysis (normally 5 mL added to the vial before shipping it to the field plus the organic-free reagent water added by the instrument). The calibration standards should also contain approximately the same amount of the sodium bisulfate preservative as the sample (e.g., ~1 g), as the presence of the preservative will affect the purging efficiencies of the analytes. The internal standard solution must be added automatically, by the instrument, in the same fashion as used for the samples. Place the soil vial containing the solution in the instrument carousel. In order to calibrate the surrogates using standards at five concentrations, it may be necessary to disable the automatic addition of surrogates to each vial containing a calibration standard (consult the manufacturer's instructions). Prior to purging, heat the sample vial to 40°C for 1.5 minutes, or as recommended by the manufacturer.

7.2.1.6 Carry out the purge-and-trap procedure as outlined in Secs. 7.2.3. to 7.2.5.

7.2.1.7 Calculate calibration factors (CF) or response factors (RF) for each analyte of interest using the procedures described in Method 8000. Calculate the average CF (external standards) or RF (internal standards) for each compound, as described in Method 8000. Evaluate the linearity of the calibration data, or choose another calibration model, as described in Method 8000 and the specific determinative method.

7.2.1.8 For GC/MS analysis, a system performance check must be made before this calibration curve is used (see Method 8260). If the purge-and-trap procedure is used with Method 8021, evaluate the response for the following four compounds: chloromethane; 1,1-dichloroethane; bromoform; and 1,1,2,2-tetrachloroethane. They are used to check for proper purge flow and to check for degradation caused by contaminated lines or active sites in the system.

7.2.1.8.1 Chloromethane is the most likely compound to be lost if the purge flow is too fast.

7.2.1.8.2 Bromoform is one of the compounds most likely to be purged very poorly if the purge flow is too slow. Cold spots and/or active sites in the transfer lines may adversely affect response.

7.2.1.8.3 Tetrachloroethane and 1,1-dichloroethane are degraded by contaminated transfer lines in purge-and-trap systems and/or active sites in trapping materials.

7.2.1.9 When analyzing for very late eluting compounds with Method 8021 (i.e., hexachlorobutadiene, 1,2,3-trichlorobenzene, etc.), cross-contamination and memory effects from a high concentration sample or even the standard are a common problem.

Extra rinsing of the purge chamber after analysis normally corrects this. The newer purge-and-trap systems often overcome this problem with better bakeout of the system following the purge-and-trap process. Also, the charcoal traps retain less moisture and decrease the problem.

7.2.2 Calibration verification

Refer to Method 8000 for details on calibration verification. A single standard near the mid-point of calibration range is used for verification. This standard should also contain approximately 1 g of sodium bisulfate.

7.2.3 Sample purge-and-trap

This method is designed for a 5-g sample size, but smaller sample sizes may be used. Consult the instrument manufacturer's instructions regarding larger sample sizes, in order to avoid clogging of the purging apparatus. The soil vial is hermetically sealed at the sampling site, and MUST remain so in order to guarantee the integrity of the sample. Gloves must be worn when handling the sample vial since the vial has been tared. If any soil is noted on the exterior of the vial or cap, it must be carefully removed prior to weighing. Weigh the vial and contents to the nearest 0.01 g, even if the sample weight was determined in the field, and record this weight. This second weighing provides a check on the field sampling procedures and provides additional assurance that the reported sample weight is accurate. Data users should be advised on significant discrepancies between the field and laboratory weights.

7.2.3.1 Remove the sample vial from storage and allow it to warm to room temperature. Shake the vial gently, to ensure that the contents move freely and that stirring will be effective. Place the sample vial in the instrument carousel according to the manufacturer's instructions.

7.2.3.2 Without disturbing the hermetic seal on the sample vial, add 5 mL of organic-free reagent water, the internal standards, and the surrogate compounds. This is carried out using the automated sampler. Other volumes of organic-free reagent water may be used, however, it is imperative that all samples, blanks, and calibration standards have exactly the same final volume of organic-free reagent water. Prior to purging, heat the sample vial to 40°C for 1.5 minutes, or as described by the manufacturer.

7.2.3.3 For the sample selected for matrix spiking, add the matrix spiking solution described in Sec. 5.0 of Method 5000, either manually, or automatically, following the manufacturer's instructions. The concentration of the spiking solution and the amount added should be established as described in Sec. 8.0 of Method 8000.

7.2.3.4 Purge the sample with helium or another inert gas at a flow rate of up to 40 mL/minute (the flow rate may vary from 20 to 40 mL/min, depending on the target analyte group) for 11 minutes while the sample is being agitated with the magnetic stirring bar or other mechanical means. The purged analytes are allowed to flow out of the vial through a glass-lined transfer line to a trap packed with suitable sorbent materials.

7.2.4 Sample Desorption

7.2.4.1 Non-cryogenic interface - After the 11 minute purge, place the purge-and-trap system in the desorb mode and preheat the trap to 245°C without a flow

of desorption gas. Start the flow of desorption gas at 10 mL/minute for about four minutes (1.5 min is normally adequate for analytes in Method 8015). Begin the temperature program of the gas chromatograph and start data acquisition.

7.2.4.2 Cryogenic interface - After the 11 minute purge, place the purge-and-trap system in the desorb mode, make sure that the cryogenic interface is at -150°C or lower, and rapidly heat the trap to 245°C while backflushing with an inert gas at 4 mL/minute for about 5 minutes (1.5 min is normally adequate for analytes in Methods 8015). At the end of the 5-minute desorption cycle, rapidly heat the cryogenic trap to 250°C. Begin the temperature program of the gas chromatograph and start the data acquisition.

7.2.5 Trap Reconditioning

After desorbing the sample for 4 minutes, recondition the trap by returning the purge-and-trap system to the purge mode. Maintain the trap temperature at 245°C (or other temperature recommended by the manufacturer of the trap packing materials). After approximately 10 minutes, turn off the trap heater and halt the purge flow through the trap. When the trap is cool, the next sample can be analyzed.

7.2.6 Data Interpretation

Perform qualitative and quantitative analysis following the guidance given in the determinative method and Method 8000. If the concentration of any target analyte exceeds the calibration range of the instrument, it will be necessary to reanalyze the sample by the high concentration method. Such reanalyses need only address those analytes for which the concentration exceeded the calibration range of the low concentration method. Alternatively, if a sample aliquot of 1-2 g was also collected (see Sec. 6.2.1.7), it may be practical to analyze that aliquot for the analytes that exceeded the instrument calibration range in the 5-g analysis. If results are to be reported on a dry weight basis, proceed to Sec. 7.5

7.3 High concentration method for soil samples with concentrations generally greater than 200 µg/kg.

The high concentration method for soil is based on a solvent extraction. A solid sample is either extracted or diluted, depending on sample solubility in a water-miscible solvent. An aliquot of the extract is added to organic-free reagent water containing surrogates and, if applicable, internal and matrix spiking standards, purged according to Method 5030, and analyzed by an appropriate determinative method. Wastes that are insoluble in methanol (i.e., petroleum and coke wastes) are diluted with hexadecane (see Sec. 7.3.8).

The specific sample preparation steps depend on whether or not the sample was preserved in the field. Samples that were not preserved in the field are prepared using the steps below, beginning at Sec. 7.3.1. If solvent preservation was employed in the field, then the preparation begins with Sec. 7.3.4.

7.3.1 When the high concentration sample is not preserved in the field, the sample consists of the entire contents of the sample container. Do not discard any supernatant liquids. Whenever practical, mix the contents of the sample container by shaking or other mechanical means without opening the vial. When shaking is not practical, quickly mix the contents of the vial with a narrow metal spatula and immediately reseal the vial.

7.3.2 If the sample is from an unknown source, perform a solubility test before proceeding. Remove several grams of material from the sample container. Quickly reseal the container to minimize the loss of volatiles. Weigh 1-g aliquots of the sample into several test tubes or other suitable containers. Add 10 mL of methanol to the first tube, 10 mL of PEG to the second, and 10 mL of hexadecane to the third. Swirl the sample and determine if it is soluble in the solvent. Once the solubility has been evaluated, discard these test solutions. If the sample is soluble in either methanol or PEG, proceed with Sec. 7.3.3. If the sample is only soluble in hexadecane, proceed with Sec. 7.3.8.

7.3.3 For soil and solid waste samples that are soluble in methanol, add 9.0 mL of methanol and 1.0 mL of the surrogate spiking solution to a tared 20-mL vial. Using a top-loading balance, weigh 5 g (wet weight) of sample into the vial. Quickly cap the vial and reweigh the vial. Record the weight to 0.1 g. Shake the vial for 2 min. If the sample was not soluble in methanol, but was soluble in PEG, employ the same procedure described above, but use 9.0 mL of PEG in place of the methanol. Proceed with Sec. 7.3.5.

NOTE: The steps in Secs. 7.3.1, 7.3.2, and 7.3.3 must be performed rapidly and without interruption to avoid loss of volatile organics. These steps must be performed in a laboratory free from solvent fumes.

7.3.4 For soil and solid waste samples that were collected in methanol or PEG (see Sec. 6.2.2), weigh the vial to 0.1 g as a check on the weight recorded in the field, add the surrogate spiking solution to the vial by injecting it through the septum, shake for 2 min, as described above, and proceed with Sec. 7.3.5.

7.3.5 Pipet approximately 1 mL of the extract from either Sec. 7.3.3 or 7.3.4 into a GC vial for storage, using a disposable pipet, and seal the vial. The remainder of the extract may be discarded. Add approximately 1 mL of methanol or PEG to a separate GC vial for use as the method blank for each set of samples extracted with the same solvent.

7.3.6 The extracts must be stored at 4°C in the dark, prior to analysis. Add an appropriate aliquot of the extract (see Table 2) to 5.0 mL of organic-free reagent water and analyze by Method 5030 in conjunction with the appropriate determinative method. Proceed to Sec. 7.0 in Method 5030 and follow the procedure for purging high concentration samples.

7.3.7 If results are to be reported on a dry weight basis, determine the dry weight of a separate aliquot of the sample, using the procedure in Sec. 7.5, after the sample extract has been transferred to a GC vial and the vial sealed.

7.3.8 For solids that are not soluble in methanol or PEG (including those samples consisting primarily of petroleum or coking waste) dilute or extract the sample with hexadecane using the procedures in Sec. 7.0 of Method 3585.

7.4 High concentration method for oily waste samples

This procedure for the analysis of oily waste samples involves the dilution of the sample in methanol or PEG. However, care must be taken to avoid introducing any of the floating oil layer into the instrument. A portion of the diluted sample is then added to 5.0 mL of organic-free reagent water, purged according to Method 5030, and analyzed using an appropriate determinative method.

For oily samples that are not soluble in methanol or PEG (including those samples consisting primarily of petroleum or coking waste), dilute or extract with hexadecane using the procedures in Sec. 7.0 of Method 3585.

The specific sample preparation steps depend on whether or not the sample was preserved in the field. Samples that were not preserved in the field are prepared using the steps below, beginning at Sec. 7.4.1. If methanol preservation was employed in the field, then the preparation begins with Sec. 7.4.3.

7.4.1 If the waste was not preserved in the field and it is soluble in methanol or PEG, weigh 1 g (wet weight) of the sample into a tared 10-mL volumetric flask, a tared scintillation vial, or a tared culture tube. If a vial or tube is used instead of a volumetric flask, it must be calibrated prior to use. This operation must be performed prior to opening the sample vial and weighing out the aliquot for analysis.

7.4.1.1 To calibrate the vessel, pipet 10.0 mL of methanol or PEG into the vial or tube and mark the bottom of the meniscus.

7.4.1.2 Discard this solvent, and proceed with weighing out the 1-g sample aliquot.

7.4.2 Quickly add 1.0 mL of surrogate spiking solution to the flask, vial, or tube, and dilute to 10.0 mL with the appropriate solvent (methanol or PEG). Swirl the vial to mix the contents and then shake vigorously for 2 minutes.

7.4.3 If the sample was collected in the field in a vial containing methanol or PEG, weigh the vial to 0.1 g as a check on the weight recorded in the field, add the surrogate spiking solution to the vial by injecting it through the septum. Swirl the vial to mix the contents and then shake vigorously for 2 minutes and proceed with Sec. 7.4.4.

7.4.4 Regardless of how the sample was collected, the target analytes are extracted into the solvent along with the majority of the oily waste (i.e., some of the oil may still be floating on the surface). If oil is floating on the surface, transfer 1 to 2 mL of the extract to a clean GC vial using a Pasteur pipet. Ensure that no oil is transferred to the vial.

7.4.5 Add 10 - 50 μ L of the methanol extract to 5 mL of organic-free reagent water for purge-and-trap analysis, using Method 5030.

7.4.6 Prepare a matrix spike sample by adding 10 - 50 μ L of the matrix spike standard dissolved in methanol to a 1-g aliquot of the oily waste. Shake the vial to disperse the matrix spike solution throughout the oil. Then add 10 mL of extraction solvent and proceed with the extraction and analysis, as described in Secs. 7.4.2 - 7.4.5. Calculate the recovery of the spiked analytes as described in Method 8000. If the recovery is not within the acceptance limits for the application, use the hexadecane dilution technique in Sec. 7.0 of Method 3585.

7.5 Determination of % Dry Weight

If results are to be reported on a dry weight basis, it is necessary to determine the dry weight of the sample.

NOTE: It is highly recommended that the dry weight determination only be made after the analyst has determined that no sample aliquots will be taken from the 60-mL vial for high

concentration analysis. This is to minimize loss of volatiles and to avoid sample contamination from the laboratory atmosphere. There is no holding time associated with the dry weight determination. Thus, this determination can be made any time prior to reporting the sample results, as long as the vial containing the additional sample has remained sealed and properly stored.

7.5.1 Weigh 5-10 g of the sample from the 60-mL VOA vial into a tared crucible.

7.5.2 Dry this aliquot overnight at 105°C. Allow to cool in a desiccator before weighing. Calculate the % dry weight as follows:

$$\% \text{ dry weight} = \frac{\text{g of dry sample}}{\text{g of sample}} \times 100$$

WARNING: The drying oven should be contained in a hood or vented. Significant laboratory contamination may result from a heavily contaminated hazardous waste sample.

8.0 QUALITY CONTROL

8.1 Refer to Chapter One for specific quality control procedures and Method 5000 for sample preparation QC procedures.

8.2 Before processing any samples, the analyst should demonstrate through the analysis of an organic-free reagent water method blank that all glassware and reagents are interference free. Each time a set of samples is extracted, or there is a change in reagents, a method blank should be processed as a safeguard against chronic laboratory contamination. The blank samples should be carried through all stages of the sample preparation and measurement.

8.3 Initial Demonstration of Proficiency - Each laboratory must demonstrate initial proficiency with each sample preparation and determinative method combination it utilizes, by generating data of acceptable accuracy and precision for target analytes in a clean matrix. The laboratory must also repeat this demonstration whenever new staff are trained or significant changes in instrumentation are made. See Sec. 8.0 of Methods 5000 and 8000 for information on how to accomplish this demonstration.

8.4 Sample Quality Control for Preparation and Analysis - See Sec. 8.0 in Method 5000 and Method 8000 for procedures to follow to demonstrate acceptable continuing performance on each set of samples to be analyzed. These include the method blank, either a matrix spike/matrix spike duplicate or a matrix spike and duplicate sample analysis, a laboratory control sample (LCS), and the addition of surrogates to each sample and QC sample.

8.5 It is recommended that the laboratory adopt additional quality assurance practices for use with this method. The specific practices that are most productive depend upon the needs of the laboratory and the nature of the samples. Whenever possible, the laboratory should analyze standard reference materials and participate in relevant performance evaluation studies.

9.0 METHOD PERFORMANCE

9.1 Single laboratory accuracy and precision data were obtained for the method analytes in three soil matrices, sand, a soil collected 10 feet below the surface of a hazardous landfill, called the

C-Horizon, and a surface garden soil. Each sample was fortified with the analytes at a concentration of 20 ng/5 g, which is equivalent to 4 µg/kg. These data are listed in tables found in Method 8260.

9.2 Single laboratory accuracy and precision data were obtained for certain method analytes when extracting oily liquid using methanol as the extraction solvent. The data are presented in a table in Method 8260. The compounds were spiked into three portions of an oily liquid (taken from a waste site) following the procedure for matrix spiking described in Sec. 7.4. This represents a worst case set of data based on recovery data from many sources of oily liquid.

10.0 REFERENCES

1. Bellar, T., "Measurement of Volatile Organic Compounds in Soils Using Modified Purge-and-Trap and Capillary Gas Chromatography/Mass Spectrometry" U.S. Environmental Protection Agency, Environmental Monitoring Systems Laboratory, Cincinnati, OH, November 1991.
2. Siegrist, R. L., Jenssen, P. D., "Evaluation of Sampling Method Effects on Volatile Organic Compound Measurements in Contaminated Soils", *Envir Sci Technol*, 1990; 24; 1387-92.
3. Hewitt, A. D., Jenkins, T. F., Grant, C. L., "Collection, Handling and Storage: Keys to Improved Data Quality for Volatile Organic Compounds in Soil", *Am Environ Lab*, 1995; 7(1); 25-8.
4. Liikala, T. L., Olsen, K. B., Teel, S. S., Lanigan, D. C., "Volatile Organic Compounds: Comparison of Two Sample Collection and Preservation Methods", *Envir Sci Technol*, 1996; 30; 3441-7.
5. Lewis, T. E., Crockett, A. B., Siegrist, R. L., Zarrabi, K., "Soil Sampling and Analysis for Volatile Organic Compounds", *Envir Monitoring & Assessment*, 1994; 30; 213-46.
6. Hewitt, A. D., "Enhanced Preservation of Volatile Organic Compounds in Soil with Sodium Bisulfate", SR95-26, U. S. Army Cold Regions Research and Engineering Laboratory, Hanover, NH.
7. Hewitt, A. D., Lukash, N. J. E., "Sampling for In-Vial Analysis of Volatile Organic Compounds in Soil", *Am Environ Lab*, 1996; Aug; 15-9.
8. Hewitt, A. D., Miyares, P. H., Sletten, R. S., "Determination of Two Chlorinated Volatile Organic Compounds in Soil by Headspace Gas Chromatography and Purge-and-Trap Gas Chromatography/Mass Spectrometry", *Hydrocarbon Contaminated Soils*, 1993, 3; 135-45, Chelsea, MI, Lewis Publishers.
9. Hewitt, A. D., "Methods of Preparing Soil Samples for Headspace Analysis of Volatile Organic Compounds: Emphasis on Salting Out", 12th Annual Waste Testing and Quality Assurance Symposium, Washington, DC, 1996, 322-9.
10. Hewitt, A. D., Miyares, P. H., Leggett, D. C., Jenkins, T. F., "Comparison of Analytical Methods for Determination of Volatile Organic Compounds", *Envir Sci Tech*, 1992; 26; 1932-8.

TABLE 1

QUANTITY OF METHANOL EXTRACT REQUIRED FOR ANALYSIS OF
HIGH CONCENTRATION SOILS/SEDIMENTS


Approximate Concentration Range				Volume of Methanol Extract ^a
500	-	10,000	µg/kg	100 µL
1,000	-	20,000	µg/kg	50 µL
5,000	-	100,000	µg/kg	10 µL
25,000	-	500,000	µg/kg	100 µL of 1/50 dilution ^b

Calculate appropriate dilution factor for concentrations exceeding those in this table.


^a The volume of methanol added to 5 mL of water being purged should be kept constant. Therefore, add to the 5-mL syringe whatever volume of methanol is necessary to maintain a total volume of 100 µL of methanol.

^b Dilute an aliquot of the methanol extract and then take 100 µL for analysis.


METHOD 5035
CLOSED-SYSTEM PURGE-AND-TRAP AND EXTRACTION
FOR VOLATILE ORGANICS IN SOIL AND WASTE SAMPLES

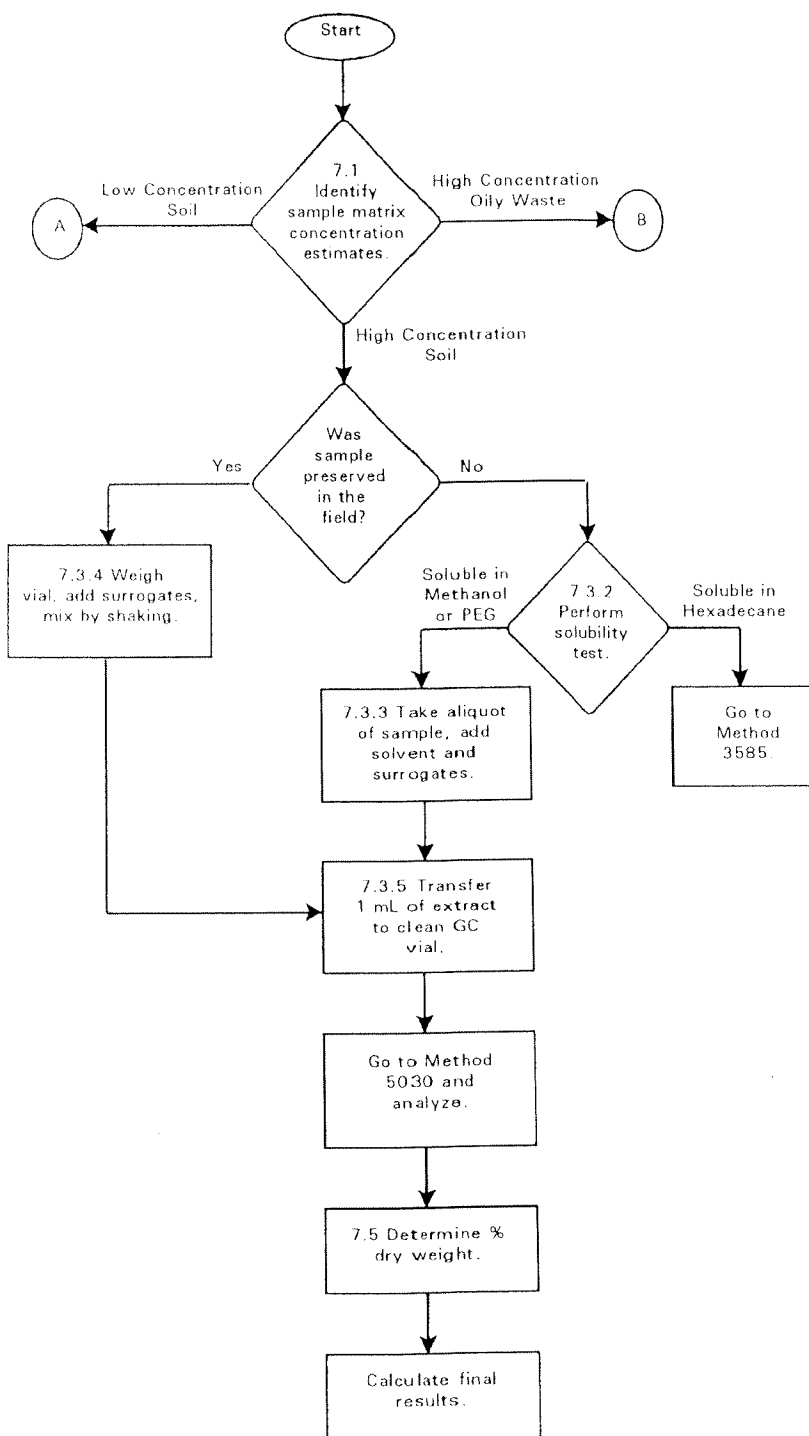

METHOD 5035 (CONTINUED)

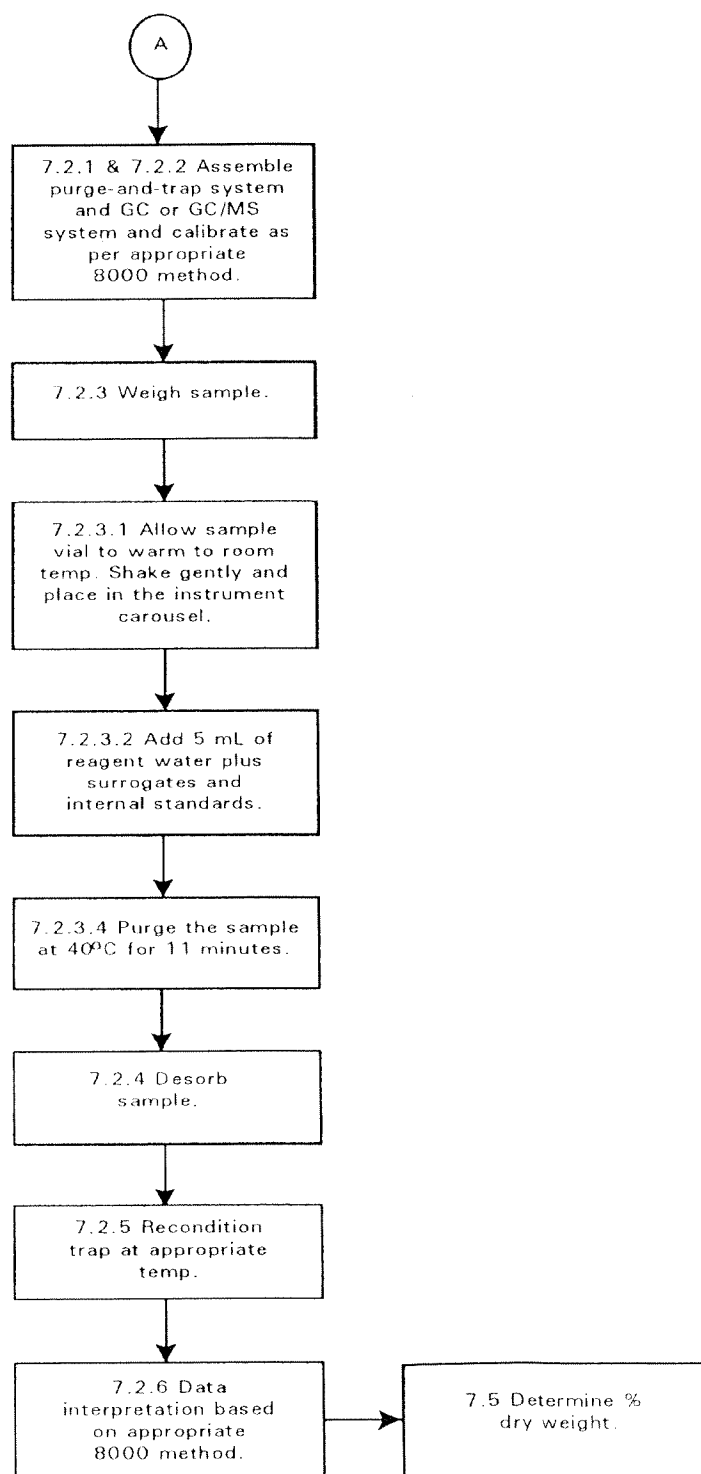

METHOD 5035 (CONTINUED)

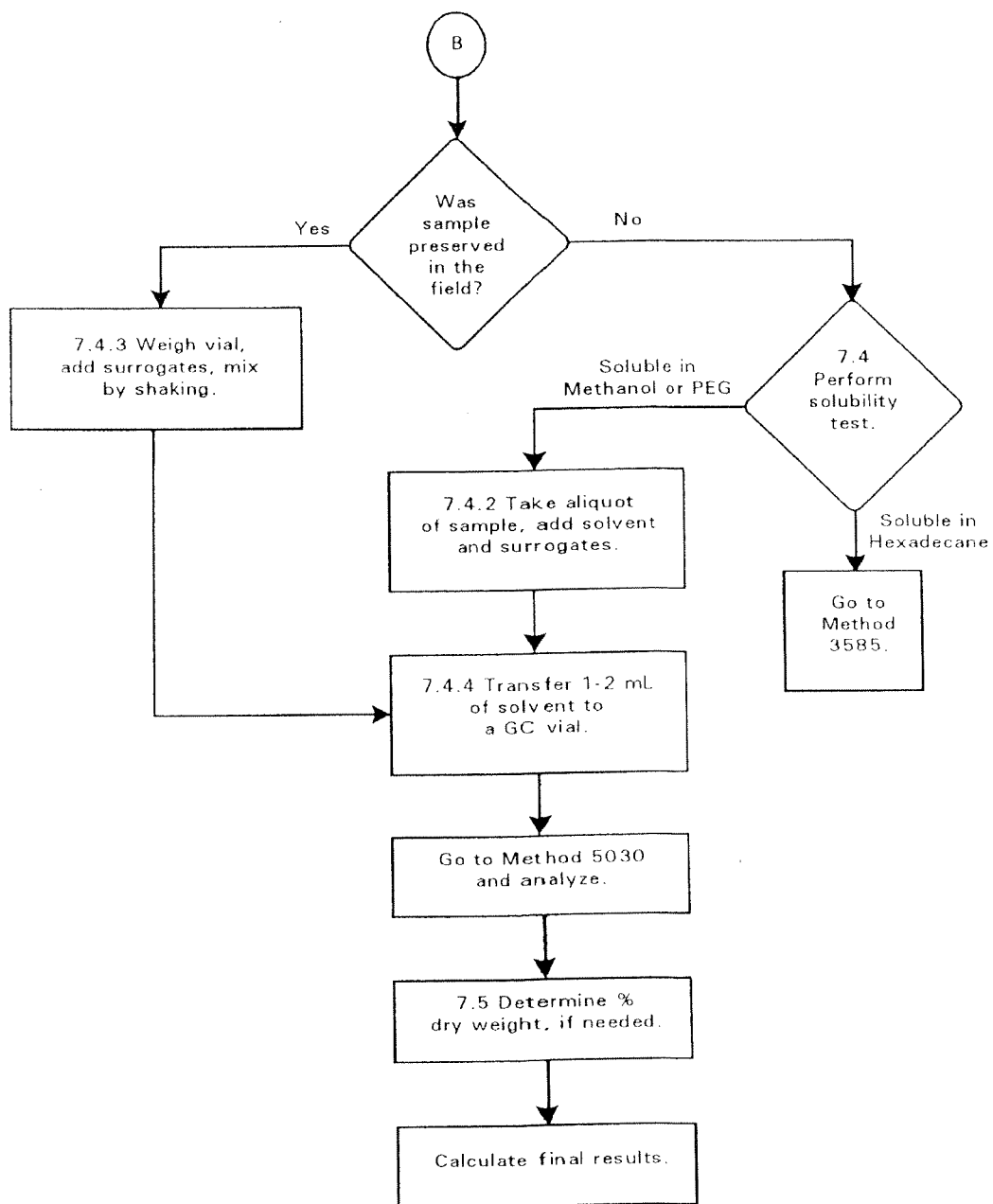

Title 22 Metals by Inductively Coupled Plasma (ICP), Graphite Furnace Atomic Absorption (GFAA), and Cold Vapor Atomic Absorption (CVAA)

SW-846 Method 6010: Inductively Coupled Plasma (ICP);
SW-846 Method 7000: Graphite Furnace Atomic Absorption (GFAA); and
SW-846 Method 7470/7471: Mercury Analysis by CVAA;

Table 1A. Summary of Holding Times and Preservation for Title 22 Metals

Analytical Parameter ^a	Technical and Contract Holding Times	Preservation
Metals in water (except Mercury)	Technical: 180 days from date of collection; Contract: 35 days from sample receipt at laboratory	pH <2 (with nitric acid)
Metals in soil (except Mercury)	Technical: 180 days from date of collection; Contract: 35 days from sample receipt at laboratory	Cool to 4°C ±2°C
Mercury in water	Technical: 28 days from date of collection; Contract: 26 days from sample receipt at laboratory	pH <2 (with nitric acid)
Mercury in soil	Technical: 28 days from date of collection; Contract: 26 days from sample receipt at laboratory	Cool to 4°C ±2°C

^a Individual target compounds are listed in Table 1B.

Sample Preparation:

Water samples are to be prepared following the protocol presented in SW-846 Method 3010. Soil samples are to be prepared following the protocol presented in SW-846 Method 3050.

Data Calculations and Reporting Units:

Calculate the sample results according to the protocol of the appropriate analytical method used.

Report water sample results in concentration units of micrograms per liter (µg/L), and soil sample results in concentrations units of milligrams per kilogram (mg/kg) on a dry weight basis. Report percent solids to the nearest percent.

For rounding results, adhere to the following rules:

- If the number following those to be retained is less than 5, round down;
- If the number following those to be retained is greater than 5, round up; or
- If the number following the last digit to be retained is equal to 5, round down if the digit is even, or round up if the digit is odd.

All records of analysis and calculations must be legible and sufficient to recalculate all sample concentrations and QC results. Include an example calculation in the data package.

TABLE 1B. Target Compound List, CAS Numbers, and Contract Required Detection Limits for Title 22 Metals by ICP, GFAA, and CVAA

COMPOUND	CAS No.	CRDL for Water (µg/L)	CRDL for Soil (mg/Kg)
Aluminum	7429-90-5	50	10.0
Antimony	7440-36-0	6	2.5
Arsenic	7440-38-2	2	1.0
Barium	7440-39-3	20	4.0
Beryllium	7440-41-7	1	0.25
Cadmium	7440-43-9	1	0.5
Chromium	7440-47-3	10	2.0
Cobalt	7440-48-4	20	4.0
Copper	7440-50-8	20	4.0
Lead	7439-92-1	5	1.0
Mercury	7439-97-6	0.2	0.1
Molybdenum	7439-97-6	20	4.0
Nickel	7440-02-0	10	4.0
Selenium	7782-49-2	5	1.0
Silver	7440-22-4	20	4.0
Thallium	7440-28-0	1	1.0
Vanadium	7440-62-2	20	4.0
Zinc	7440-66-6	20	5.0

Table 2A. Summary of Calibration Procedures for Title 22 Metals by ICP

Calibration Element	Frequency	Acceptance Criteria	Corrective Action
Initial Calibration (minimum blank + 1 calibration standard) (ICAL)	Initially, Daily; whenever required, due to failure of CCV	Acceptable ICV, CRDL, and ICB standards	1. Terminate analysis 2. Re-calibrate and verify before sample analysis
Initial Calibration Verification (ICV) at midpoint of ICAL (Different source from ICAL standards)	Daily, immediately following ICAL and prior to sample analysis	$\pm 10\%$ from expected concentration	1. Terminate analysis and identify and document problem 2. Reprep and re-analyze ICV and all associated samples 3. Re-calibrate and re-analyze reprepped ICV and all associated samples
Calibration Blank Verification (ICB, CCB)	After ICV and every CCV	< CRDL	1. Terminate analysis 2. Determine Source of contamination 3. Reprep ICB and CCB 4. Re-analyze all samples associated with a contaminated blank
Continuing Calibration Verification (CCV)	Before samples, after every 10 samples, and end of run	$\pm 10\%$ from expected concentration	1. Re-calibrate and verify 2. Re-analyze samples back to last acceptable CCV
Contract Required Detection Limit Verification Standard (CRI)	After ICV, but before sample analysis	$\pm 35\%$ from expected concentration	1. Re-calibrate and verify 2. Re-analyze samples back to last compliant CCV
ICP Interference Check Sample (ICS)	Run at start and finish of daily run or twice per 8 hours	$\pm 20\%$ from true value concentration	1. Reprep and re-analyze standard 2. Re-calibrate, verify and re-analyze all associated samples

Table 2B. Summary of Calibration Procedures for Title 22 Metals by GFAA

Calibration Element	Frequency	Acceptance Criteria	Corrective Action
Initial Calibration (minimum blank + 3 standards for each analyte) (ICAL) ^a	Initially, each analytical batch; whenever required, due to failure of CCV	$r \geq 0.995$	1. Terminate analysis 2. Re-calibrate and verify before sample analysis
Initial Calibration Verification (ICV) at midpoint of ICAL (Different source from ICAL standards)	Daily, immediately following ICAL and prior to sample analysis	$\pm 10\%$ from expected concentration	1. Terminate analysis and identify and document problem 2. Reprep and re-analyze ICV and all associated samples 3. Re-calibrate and re-analyze reprepped ICV and all associated samples
Calibration Blank Verification (ICB, CCB)	After ICV and every CCV	$< \text{CRDL}$	1. Terminate analysis 2. Determine Source of contamination 3. Reprep ICB and CCB 4. Re-analyze all samples associated with a contaminated blank
Continuing Calibration Verification (CCV)	Following ICV and before sample analysis; after every 10 samples and end of run	$\pm 10\%$ from expected concentration	1. Re-calibrate and verify 2. Re-analyze samples back to last acceptable CCV
Contract Required Detection Limit Verification Standard (CRA)	After ICV, but before sample analysis	$\pm 35\%$ from expected concentration	1. Reprep and re-analyze standard 2. Re-calibrate and verify

^a The ICAL low standard must be between the CRDL and 2X CRDL.

Table 2C. Summary of Calibration Procedures for Mercury by CVAA

Calibration Element	Frequency	Acceptance Criteria	Corrective Action
Initial Calibration (minimum blank + 5 standards) (ICAL) ^a	Initially, each analytical batch; whenever required, due to failure of CCV	$r \geq 0.995$	1. Terminate analysis 2. Re-calibrate and verify before sample analysis
Initial Calibration Verification (ICV) at midpoint of ICAL (Different source from ICAL standards)	Daily, immediately following ICAL and prior to sample analysis	$\pm 20\%$ from expected concentration	1. Terminate analysis and identify and document problem 2. Reprep and re-analyze ICV and all associated samples 3. Re-calibrate and re-analyze reprepped ICV and all associated samples
Calibration Blank Verification (ICB, CCB)	After ICV and every CCV	$< \text{CRDL}$	1. Terminate analysis 2. Determine Source of contamination 3. Reprep ICB and CCB 4. Re-analyze all samples associated with a contaminated blank
Continuing Calibration Verification (CCV)	Before Samples, after every 10 samples, and end of run	$\pm 20\%$ from expected concentration	1. Re-calibrate and verify 2. Re-analyze samples back to last acceptable CCV
Contract Required Detection Limit Verification Standard (CRA)	After ICV, but before sample analysis	$\pm 35\%$ from expected concentration	1. Reprep and re-analyze standard 2. Re-calibrate and verify

^a The ICAL low standard must be at the CRDL.

Table 3. Summary of Internal Quality Control Procedures for Title 22 Metals by ICP, CVAA, and GFAA

QC Element	Frequency	Acceptance Criteria	Corrective Action
Method Blank (MB)	One per batch or SDG ^{a, b}	< CRDL	1. If lowest sample concentration is more than 10X the blank conc., no action 2. If samples are non-detected, no action 3. If detected sample concentrations are less than 10X blank conc., all affected samples must be prepared again with another method blank and re-analyzed
Duplicate Sample (DUP)	One per batch or SDG ^{a, b}	Waters: RPD <± 20% for samples >5X CRDL; ± CRDL for samples <5X CRDL Soils: RPD <± 35% for samples >5X CRDL; ± 2xCRDL for samples <5X CRDL	1. Flag associated data with an "*"
Matrix Spike Sample (MS)	One per batch or SDG ^{a, b}	± 25% from expected value ^c	1. Flag associated data with an "N"
Laboratory Control Sample (LCS)	One per batch or SDG ^{a, b}	Waters: ± 20% from expected concentration Soils: within control limits of certified solid LCS or ± 20% from expected spike concentration	1. Terminate analysis and identify and document the problem 2. Re-analyze all associated samples
Serial Dilution Sample (5 X Dilution) (ICP only)	One per batch or SDG ^{a, b}	± 10% difference from original results for analytes greater than 50 X IDL	1. Flag associated data with a "B"

Table 3. (cont) Summary of Internal Quality Control Procedures for Title 22 Metals by ICP, CVAA, and GFAA

QC Element	Frequency	Acceptance Criteria	Corrective Action
Duplicate Injections (GFAA only)	All samples	Duplicate results within $\pm 20\%$ RPD (or CV)	1. Rerun sample once 2. Flag associated data with an "E" if acceptance criteria are not met after second run
Analytical Spike Sample (2 X CRDL) (GFAA only)	All samples	Spike recovery $\pm 15\%$	1. If spike recovery is $<40\%$, dilute sample by a factor of 5 to 10 and run again. If recovery is still $<40\%$, report data and flag with an "E" to indicate interference problems 2. If sample concentration is $<50\%$ of recovered spike value and spike recovery is $>40\%$ and $<85\%$ or $>115\%$, report result down to IDL and flag result with a "W" 3. If sample concentration is $>50\%$ of recovered spike value and spike recovery is $>40\%$ and $<85\%$ or $>115\%$, quantitate by MSA
Method of Standard Addition (GFAA only)	As determined by analytical spike recovery results	3 samples spiked at 50%, 100%, and 150% of sample concentration and $r \geq 0.995$	1. Rerun samples only once 2. Flag associated data with a "+" if acceptance criteria are not met after second run

^a SDG - Sample Delivery Group - each case of field samples received; or each 20 field samples within a case; or each 14 calendar day period during which field samples in a case are received.

^b Minimum requirement is the analysis of 1 QC sample per 20 samples.

^c An exception to this rule is granted in situations where the sample concentration exceeds the spike concentration by a factor of 4. In such an event, the data shall be reported unflagged.

Dilute and reanalyze samples with concentrations exceeding the range of the calibration curve. Results for such re-analyses should fall within the mid-range of the calibration curve. Report results and submit documentation for both analyses.

METHOD 6010B

INDUCTIVELY COUPLED PLASMA-ATOMIC EMISSION SPECTROMETRY

1.0 SCOPE AND APPLICATION

1.1 Inductively coupled plasma-atomic emission spectrometry (ICP-AES) determines trace elements, including metals, in solution. The method is applicable to all of the elements listed in Table 1. All matrices, excluding filtered groundwater samples but including ground water, aqueous samples, TCLP and EP extracts, industrial and organic wastes, soils, sludges, sediments, and other solid wastes, require digestion prior to analysis. Groundwater samples that have been prefiltered and acidified will not need acid digestion. Samples which are not digested must either use an internal standard or be matrix matched with the standards. Refer to Chapter Three for the appropriate digestion procedures.

1.2 Table 1 lists the elements for which this method is applicable. Detection limits, sensitivity, and the optimum and linear concentration ranges of the elements can vary with the wavelength, spectrometer, matrix and operating conditions. Table 1 lists the recommended analytical wavelengths and estimated instrumental detection limits for the elements in clean aqueous matrices. The instrument detection limit data may be used to estimate instrument and method performance for other sample matrices. Elements and matrices other than those listed in Table 1 may be analyzed by this method if performance at the concentration levels of interest (see Section 8.0) is demonstrated.

1.3 Users of the method should state the data quality objectives prior to analysis and must document and have on file the required initial demonstration performance data described in the following sections prior to using the method for analysis.

1.4 Use of this method is restricted to spectroscopists who are knowledgeable in the correction of spectral, chemical, and physical interferences described in this method.

2.0 SUMMARY OF METHOD

2.1 Prior to analysis, samples must be solubilized or digested using appropriate Sample Preparation Methods (e.g. Chapter Three). When analyzing groundwater samples for dissolved constituents, acid digestion is not necessary if the samples are filtered and acid preserved prior to analysis.

2.2 This method describes multielemental determinations by ICP-AES using sequential or simultaneous optical systems and axial or radial viewing of the plasma. The instrument measures characteristic emission spectra by optical spectrometry. Samples are nebulized and the resulting aerosol is transported to the plasma torch. Element-specific emission spectra are produced by a radio-frequency inductively coupled plasma. The spectra are dispersed by a grating spectrometer, and the intensities of the emission lines are monitored by photosensitive devices. Background correction is required for trace element determination. Background must be measured adjacent to analyte lines on samples during analysis. The position selected for the background-intensity measurement, on either or both sides of the analytical line, will be determined by the complexity of the spectrum adjacent to the analyte line. In one mode of analysis the position used should be as free as possible from spectral interference and should reflect the same change in background

intensity as occurs at the analyte wavelength measured. Background correction is not required in cases of line broadening where a background correction measurement would actually degrade the analytical result. The possibility of additional interferences named in Section 3.0 should also be recognized and appropriate corrections made; tests for their presence are described in Section 8.5. Alternatively, users may choose multivariate calibration methods. In this case, point selections for background correction are superfluous since whole spectral regions are processed.

3.0 INTERFERENCES

3.1 Spectral interferences are caused by background emission from continuous or recombination phenomena, stray light from the line emission of high concentration elements, overlap of a spectral line from another element, or unresolved overlap of molecular band spectra.

3.1.1 Background emission and stray light can usually be compensated for by subtracting the background emission determined by measurements adjacent to the analyte wavelength peak. Spectral scans of samples or single element solutions in the analyte regions may indicate when alternate wavelengths are desirable because of severe spectral interference. These scans will also show whether the most appropriate estimate of the background emission is provided by an interpolation from measurements on both sides of the wavelength peak or by measured emission on only one side. The locations selected for the measurement of background intensity will be determined by the complexity of the spectrum adjacent to the wavelength peak. The locations used for routine measurement must be free of off-line spectral interference (interelement or molecular) or adequately corrected to reflect the same change in background intensity as occurs at the wavelength peak. For multivariate methods using whole spectral regions, background scans should be included in the correction algorithm. Off-line spectral interferences are handled by including spectra on interfering species in the algorithm.

3.1.2 To determine the appropriate location for off-line background correction, the user must scan the area on either side adjacent to the wavelength and record the apparent emission intensity from all other method analytes. This spectral information must be documented and kept on file. The location selected for background correction must be either free of off-line interelement spectral interference or a computer routine must be used for automatic correction on all determinations. If a wavelength other than the recommended wavelength is used, the analyst must determine and document both the overlapping and nearby spectral interference effects from all method analytes and common elements and provide for their automatic correction on all analyses. Tests to determine spectral interference must be done using analyte concentrations that will adequately describe the interference. Normally, 100 mg/L single element solutions are sufficient; however, for analytes such as iron that may be found at high concentration, a more appropriate test would be to use a concentration near the upper analytical range limit.

3.1.3 Spectral overlaps may be avoided by using an alternate wavelength or can be compensated by equations that correct for interelement contributions. Instruments that use equations for interelement correction **require** the interfering elements be analyzed at the same time as the element of interest. When operative and uncorrected, interferences will produce false positive determinations and be reported as analyte concentrations. More extensive information on interferant effects at various wavelengths and resolutions is available in reference wavelength tables and books. Users may apply interelement

correction equations determined on their instruments with tested concentration ranges to compensate (off line or on line) for the effects of interfering elements. Some potential spectral interferences observed for the recommended wavelengths are given in Table 2. For multivariate methods using whole spectral regions, spectral interferences are handled by including spectra of the interfering elements in the algorithm. The interferences listed are only those that occur between method analytes. Only interferences of a direct overlap nature are listed. These overlaps were observed with a single instrument having a working resolution of 0.035 nm.

3.1.4 When using interelement correction equations, the interference may be expressed as analyte concentration equivalents (i.e. false analyte concentrations) arising from 100 mg/L of the interference element. For example, assume that As is to be determined (at 193.696 nm) in a sample containing approximately 10 mg/L of Al. According to Table 2, 100 mg/L of Al would yield a false signal for As equivalent to approximately 1.3 mg/L. Therefore, the presence of 10 mg/L of Al would result in a false signal for As equivalent to approximately 0.13 mg/L. The user is cautioned that other instruments may exhibit somewhat different levels of interference than those shown in Table 2. The interference effects must be evaluated for each individual instrument since the intensities will vary.

3.1.5 Interelement corrections will vary for the same emission line among instruments because of differences in resolution, as determined by the grating, the entrance and exit slit widths, and by the order of dispersion. Interelement corrections will also vary depending upon the choice of background correction points. Selecting a background correction point where an interfering emission line may appear should be avoided when practical. Interelement corrections that constitute a major portion of an emission signal may not yield accurate data. Users should not forget that some samples may contain uncommon elements that could contribute spectral interferences.

3.1.6 The interference effects must be evaluated for each individual instrument whether configured as a sequential or simultaneous instrument. For each instrument, intensities will vary not only with optical resolution but also with operating conditions (such as power, viewing height and argon flow rate). When using the recommended wavelengths, the analyst is required to determine and document for each wavelength the effect from referenced interferences (Table 2) as well as any other suspected interferences that may be specific to the instrument or matrix. The analyst is encouraged to utilize a computer routine for automatic correction on all analyses.

3.1.7 Users of sequential instruments must verify the absence of spectral interference by scanning over a range of 0.5 nm centered on the wavelength of interest for several samples. The range for lead, for example, would be from 220.6 to 220.1 nm. This procedure must be repeated whenever a new matrix is to be analyzed and when a new calibration curve using different instrumental conditions is to be prepared. Samples that show an elevated background emission across the range may be background corrected by applying a correction factor equal to the emission adjacent to the line or at two points on either side of the line and interpolating between them. An alternate wavelength that does not exhibit a background shift or spectral overlap may also be used.

3.1.8 If the correction routine is operating properly, the determined apparent analyte(s) concentration from analysis of each interference solution should fall within a specific concentration range around the calibration blank. The concentration range is calculated by multiplying the concentration of the interfering element by the value of the correction factor being tested and divided by 10. If after the subtraction of the calibration blank the apparent analyte concentration falls outside of this range in either a positive or negative direction, a change in the correction factor of more than 10% should be suspected. The cause of the change should be determined and corrected and the correction factor updated. The interference check solutions should be analyzed more than once to confirm a change has occurred. Adequate rinse time between solutions and before analysis of the calibration blank will assist in the confirmation.

3.1.9 When interelement corrections are applied, their accuracy should be verified, daily, by analyzing spectral interference check solutions. If the correction factors or multivariate correction matrices tested on a daily basis are found to be within the 20% criteria for 5 consecutive days, the required verification frequency of those factors in compliance may be extended to a weekly basis. Also, if the nature of the samples analyzed is such they do not contain concentrations of the interfering elements at \pm one reporting limit from zero, daily verification is not required. All interelement spectral correction factors or multivariate correction matrices must be verified and updated every six months or when an instrumentation change, such as in the torch, nebulizer, injector, or plasma conditions occurs. Standard solution should be inspected to ensure that there is no contamination that may be perceived as a spectral interference.

3.1.10 When interelement corrections are not used, verification of absence of interferences is required.

3.1.10.1 One method is to use a computer software routine for comparing the determinative data to limits files for notifying the analyst when an interfering element is detected in the sample at a concentration that will produce either an apparent false positive concentration, (i.e., greater than) the analyte instrument detection limit, or false negative analyte concentration, (i.e., less than the lower control limit of the calibration blank defined for a 99% confidence interval).

3.1.10.2 Another method is to analyze an Interference Check Solution(s) which contains similar concentrations of the major components of the samples (>10 mg/L) on a continuing basis to verify the absence of effects at the wavelengths selected. These data must be kept on file with the sample analysis data. If the check solution confirms an operative interference that is \geq 20% of the analyte concentration, the analyte must be determined using (1) analytical and background correction wavelengths (or spectral regions) free of the interference, (2) by an alternative wavelength, or (3) by another documented test procedure.

3.2 Physical interferences are effects associated with the sample nebulization and transport processes. Changes in viscosity and surface tension can cause significant inaccuracies, especially in samples containing high dissolved solids or high acid concentrations. If physical interferences are present, they must be reduced by diluting the sample or by using a peristaltic pump, by using an internal standard or by using a high solids nebulizer. Another problem that can occur with high dissolved solids is salt buildup at the tip of the nebulizer, affecting aerosol flow rate

and causing instrumental drift. The problem can be controlled by wetting the argon prior to nebulization, using a tip washer, using a high solids nebulizer or diluting the sample. Also, it has been reported that better control of the argon flow rate, especially to the nebulizer, improves instrument performance: this may be accomplished with the use of mass flow controllers. The test described in Section 8.5.1 will help determine if a physical interference is present.

3.3 Chemical interferences include molecular compound formation, ionization effects, and solute vaporization effects. Normally, these effects are not significant with the ICP technique, but if observed, can be minimized by careful selection of operating conditions (incident power, observation position, and so forth), by buffering of the sample, by matrix matching, and by standard addition procedures. Chemical interferences are highly dependent on matrix type and the specific analyte element.

3.4 Memory interferences result when analytes in a previous sample contribute to the signals measured in a new sample. Memory effects can result from sample deposition on the uptake tubing to the nebulizer and from the build up of sample material in the plasma torch and spray chamber. The site where these effects occur is dependent on the element and can be minimized by flushing the system with a rinse blank between samples. The possibility of memory interferences should be recognized within an analytical run and suitable rinse times should be used to reduce them. The rinse times necessary for a particular element must be estimated prior to analysis. This may be achieved by aspirating a standard containing elements at a concentration ten times the usual amount or at the top of the linear dynamic range. The aspiration time for this sample should be the same as a normal sample analysis period, followed by analysis of the rinse blank at designated intervals. The length of time required to reduce analyte signals to within a factor of two of the method detection limit should be noted. Until the required rinse time is established, this method suggests a rinse period of at least 60 seconds between samples and standards. If a memory interference is suspected, the sample must be reanalyzed after a rinse period of sufficient length. Alternate rinse times may be established by the analyst based upon their DQOs.

3.5 Users are advised that high salt concentrations can cause analyte signal suppressions and confuse interference tests. If the instrument does not display negative values, fortify the interference check solution with the elements of interest at 0.5 to 1 mg/L and measure the added standard concentration accordingly. Concentrations should be within 20% of the true spiked concentration or dilution of the samples will be necessary. In the absence of measurable analyte, overcorrection could go undetected if a negative value is reported as zero.

3.6 The dashes in Table 2 indicate that no measurable interferences were observed even at higher interferant concentrations. Generally, interferences were discernible if they produced peaks, or background shifts, corresponding to 2 to 5% of the peaks generated by the analyte concentrations.

4.0 APPARATUS AND MATERIALS

4.1 Inductively coupled argon plasma emission spectrometer:

4.1.1 Computer-controlled emission spectrometer with background correction.

4.1.2 Radio-frequency generator compliant with FCC regulations.

- 4.1.3 Optional mass flow controller for argon nebulizer gas supply.
- 4.1.4 Optional peristaltic pump.
- 4.1.5 Optional Autosampler.
- 4.1.6 Argon gas supply - high purity.
- 4.2 Volumetric flasks of suitable precision and accuracy.
- 4.3 Volumetric pipets of suitable precision and accuracy.

5.0 REAGENTS

5.1 Reagent or trace metals grade chemicals shall be used in all tests. Unless otherwise indicated, it is intended that all reagents shall conform to the specifications of the Committee on Analytical Reagents of the American Chemical Society, where such specifications are available. Other grades may be used, provided it is first ascertained that the reagent is of sufficiently high purity to permit its use without lessening the accuracy of the determination. If the purity of a reagent is in question analyze for contamination. If the concentration of the contamination is less than the MDL then the reagent is acceptable.

5.1.1 Hydrochloric acid (conc), HCl.

5.1.2 Hydrochloric acid (1:1), HCl. Add 500 mL concentrated HCl to 400 mL water and dilute to 1 liter in an appropriately sized beaker.

5.1.3 Nitric acid (conc), HNO_3 .

5.1.4 Nitric acid (1:1), HNO_3 . Add 500 mL concentrated HNO_3 to 400 mL water and dilute to 1 liter in an appropriately sized beaker.

5.2 Reagent Water. All references to water in the method refer to reagent water unless otherwise specified. Reagent water will be interference free. Refer to Chapter One for a definition of reagent water.

5.3 Standard stock solutions may be purchased or prepared from ultra- high purity grade chemicals or metals (99.99% pure or greater). All salts must be dried for 1 hour at 105°C , unless otherwise specified.

Note: This section does not apply when analyzing samples that have been prepared by Method 3040.

CAUTION: Many metal salts are extremely toxic if inhaled or swallowed. Wash hands thoroughly after handling.

Typical stock solution preparation procedures follow. Concentrations are calculated based upon the weight of pure metal added, or with the use of the element fraction and the weight of the metal salt added.

For metals:

$$\text{Concentration (ppm)} = \frac{\text{weight (mg)}}{\text{volume (L)}}$$

For metal salts:

$$\text{Concentration (ppm)} = \frac{\text{weight (mg)} \times \text{mole fraction}}{\text{volume (L)}}$$

5.3.1 Aluminum solution, stock, 1 mL = 1000 µg Al: Dissolve 1.000 g of aluminum metal, weighed accurately to at least four significant figures, in an acid mixture of 4.0 mL of (1:1) HCl and 1.0 mL of concentrated HNO₃ in a beaker. Warm beaker slowly to effect solution. When dissolution is complete, transfer solution quantitatively to a 1-liter flask, add an additional 10.0 mL of (1:1) HCl and dilute to volume with reagent water.

NOTE: Weight of analyte is expressed to four significant figures for consistency with the weights below because rounding to two decimal places can contribute up to 4 % error for some of the compounds.

5.3.2 Antimony solution, stock, 1 mL = 1000 µg Sb: Dissolve 2.6673 g K(SbO)C₄H₄O₆ (element fraction Sb = 0.3749), weighed accurately to at least four significant figures, in water, add 10 mL (1:1) HCl, and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.3 Arsenic solution, stock, 1 mL = 1000 µg As: Dissolve 1.3203 g of As₂O₃ (element fraction As = 0.7574), weighed accurately to at least four significant figures, in 100 mL of water containing 0.4 g NaOH. Acidify the solution with 2 mL concentrated HNO₃ and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.4 Barium solution, stock, 1 mL = 1000 µg Ba: Dissolve 1.5163 g BaCl₂ (element fraction Ba = 0.6595), dried at 250°C for 2 hours, weighed accurately to at least four significant figures, in 10 mL water with 1 mL (1:1) HCl. Add 10.0 mL (1:1) HCl and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.5 Beryllium solution, stock, 1 mL = 1000 µg Be: Do not dry. Dissolve 19.6463 g BeSO₄·4H₂O (element fraction Be = 0.0509), weighed accurately to at least four significant figures, in water, add 10.0 mL concentrated HNO₃, and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.6 Boron solution, stock, 1 mL = 1000 µg B: Do not dry. Dissolve 5.716 g anhydrous H₃BO₃ (B fraction = 0.1749), weighed accurately to at least four significant figures, in reagent water and dilute in a 1-L volumetric flask with reagent water. Transfer immediately after mixing in a clean polytetrafluoroethylene (PTFE) bottle to minimize any leaching of boron from the glass volumetric container. Use of a non-glass volumetric flask is recommended to avoid boron contamination from glassware.

5.3.7 Cadmium solution, stock, 1 mL = 1000 µg Cd: Dissolve 1.1423 g CdO (element fraction Cd = 0.8754), weighed accurately to at least four significant figures, in a

minimum amount of (1:1) HNO_3 . Heat to increase rate of dissolution. Add 10.0 mL concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.8 Calcium solution, stock, 1 mL = 1000 μg Ca: Suspend 2.4969 g CaCO_3 (element Ca fraction = 0.4005), dried at 180°C for 1 hour before weighing, weighed accurately to at least four significant figures, in water and dissolve cautiously with a minimum amount of (1:1) HNO_3 . Add 10.0 mL concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.9 Chromium solution, stock, 1 mL = 1000 μg Cr: Dissolve 1.9231 g CrO_3 (element fraction Cr = 0.5200), weighed accurately to at least four significant figures, in water. When solution is complete, acidify with 10 mL concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.10 Cobalt solution, stock, 1 mL = 1000 μg Co: Dissolve 1.00 g of cobalt metal, weighed accurately to at least four significant figures, in a minimum amount of (1:1) HNO_3 . Add 10.0 mL (1:1) HCl and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.11 Copper solution, stock, 1 mL = 1000 μg Cu: Dissolve 1.2564 g CuO (element fraction Cu = 0.7989), weighed accurately to at least four significant figures, in a minimum amount of (1:1) HNO_3 . Add 10.0 mL concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.12 Iron solution, stock, 1 mL = 1000 μg Fe: Dissolve 1.4298 g Fe_2O_3 (element fraction Fe = 0.6994), weighed accurately to at least four significant figures, in a warm mixture of 20 mL (1:1) HCl and 2 mL of concentrated HNO_3 . Cool, add an additional 5.0 mL of concentrated HNO_3 , and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.13 Lead solution, stock, 1 mL = 1000 μg Pb: Dissolve 1.5985 g $\text{Pb}(\text{NO}_3)_2$ (element fraction Pb = 0.6256), weighed accurately to at least four significant figures, in a minimum amount of (1:1) HNO_3 . Add 10 mL (1:1) HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.14 Lithium solution, stock, 1 mL = 1000 μg Li: Dissolve 5.3248 g lithium carbonate (element fraction Li = 0.1878), weighed accurately to at least four significant figures, in a minimum amount of (1:1) HCl and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.15 Magnesium solution, stock, 1 mL = 1000 μg Mg: Dissolve 1.6584 g MgO (element fraction Mg = 0.6030), weighed accurately to at least four significant figures, in a minimum amount of (1:1) HNO_3 . Add 10.0 mL (1:1) concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.16 Manganese solution, stock, 1 mL = 1000 μg Mn: Dissolve 1.00 g of manganese metal, weighed accurately to at least four significant figures, in acid mixture (10 mL concentrated HCl and 1 mL concentrated HNO_3) and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.17 Mercury solution, stock, 1 mL = 1000 µg Hg: Do not dry, highly toxic element. Dissolve 1.354 g HgCl_2 (Hg fraction = 0.7388) in reagent water. Add 50.0 mL concentrated HNO_3 and dilute to volume in 1-L volumetric flask with reagent water.

5.3.18 Molybdenum solution, stock, 1 mL = 1000 µg Mo: Dissolve 1.7325 g $(\text{NH}_4)_6\text{Mo}_7\text{O}_{24} \cdot 4\text{H}_2\text{O}$ (element fraction Mo = 0.5772), weighed accurately to at least four significant figures, in water and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.19 Nickel solution, stock, 1 mL = 1000 µg Ni: Dissolve 1.00 g of nickel metal, weighed accurately to at least four significant figures, in 10.0 mL hot concentrated HNO_3 , cool, and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.20 Phosphate solution, stock, 1 mL = 1000 µg P: Dissolve 4.3937 g anhydrous KH_2PO_4 (element fraction P = 0.2276), weighed accurately to at least four significant figures, in water. Dilute to volume in a 1,000 mL volumetric flask with water.

5.3.21 Potassium solution, stock, 1 mL = 1000 µg K: Dissolve 1.9069 g KCl (element fraction K = 0.5244) dried at 110°C, weighed accurately to at least four significant figures, in water, and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.22 Selenium solution, stock, 1 mL = 1000 µg Se: Do not dry. Dissolve 1.6332 g H_2SeO_3 (element fraction Se = 0.6123), weighed accurately to at least four significant figures, in water and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.23 Silica solution, stock, 1 mL = 1000 µg SiO_2 : Do not dry. Dissolve 2.964 g NH_4SiF_6 , weighed accurately to at least four significant figures, in 200 mL (1:20) HCl with heating at 85°C to effect dissolution. Let solution cool and dilute to volume in a 1-L volumetric flask with reagent water.

5.3.24 Silver solution, stock, 1 mL = 1000 µg Ag: Dissolve 1.5748 g AgNO_3 (element fraction Ag = 0.6350), weighed accurately to at least four significant figures, in water and 10 mL concentrated HNO_3 . Dilute to volume in a 1,000 mL volumetric flask with water.

5.3.25 Sodium solution, stock, 1 mL = 1000 µg Na: Dissolve 2.5419 g NaCl (element fraction Na = 0.3934), weighed accurately to at least four significant figures, in water. Add 10.0 mL concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.26 Strontium solution, stock, 1 mL = 1000 µg Sr: Dissolve 2.4154 g of strontium nitrate ($\text{Sr}(\text{NO}_3)_2$) (element fraction Sr = 0.4140), weighed accurately to at least four significant figures, in a 1-liter flask containing 10 mL of concentrated HCl and 700 mL of water. Dilute to volume in a 1,000 mL volumetric flask with water.

5.3.27 Thallium solution, stock, 1 mL = 1000 µg Tl: Dissolve 1.3034 g TlNO_3 (element fraction Tl = 0.7672), weighed accurately to at least four significant figures, in water. Add 10.0 mL concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.28 Tin solution, stock, 1 mL = 1000 µg Sn: Dissolve 1.000 g Sn shot, weighed accurately to at least 4 significant figures, in 200 mL (1:1) HCl with heating to effect dissolution. Let solution cool and dilute with (1:1) HCl in a 1-L volumetric flask.

5.3.29 Vanadium solution, stock, 1 mL = 1000 µg V: Dissolve 2.2957 g NH_4VO_3 (element fraction V = 0.4356), weighed accurately to at least four significant figures, in a minimum amount of concentrated HNO_3 . Heat to increase rate of dissolution. Add 10.0 mL concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.3.30 Zinc solution, stock, 1 mL = 1000 µg Zn: Dissolve 1.2447 g ZnO (element fraction Zn = 0.8034), weighed accurately to at least four significant figures, in a minimum amount of dilute HNO_3 . Add 10.0 mL concentrated HNO_3 and dilute to volume in a 1,000 mL volumetric flask with water.

5.4 Mixed calibration standard solutions - Prepare mixed calibration standard solutions by combining appropriate volumes of the stock solutions in volumetric flasks (see Table 3). Add the appropriate types and volumes of acids so that the standards are matrix matched with the sample digestates. Prior to preparing the mixed standards, each stock solution should be analyzed separately to determine possible spectral interference or the presence of impurities. Care should be taken when preparing the mixed standards to ensure that the elements are compatible and stable together. Transfer the mixed standard solutions to FEP fluorocarbon or previously unused polyethylene or polypropylene bottles for storage. Fresh mixed standards should be prepared, as needed, with the realization that concentration can change on aging. Some typical calibration standard combinations are listed in Table 3.

NOTE: If the addition of silver to the recommended acid combination results in an initial precipitation, add 15 mL of water and warm the flask until the solution clears. Cool and dilute to 100 mL with water. For this acid combination, the silver concentration should be limited to 2 mg/L. Silver under these conditions is stable in a tap-water matrix for 30 days. Higher concentrations of silver require additional HCl.

5.5 Two types of blanks are required for the analysis for samples prepared by any method other than 3040. The calibration blank is used in establishing the analytical curve, and the method blank is used to identify possible contamination resulting from varying amounts of the acids used in the sample processing.

5.5.1 The calibration blank is prepared by acidifying reagent water to the same concentrations of the acids found in the standards and samples. Prepare a sufficient quantity to flush the system between standards and samples. The calibration blank will also be used for all initial and continuing calibration blank determinations (see Sections 7.3 and 7.4).

5.5.2 The method blank must contain all of the reagents in the same volumes as used in the processing of the samples. The method blank must be carried through the complete procedure and contain the same acid concentration in the final solution as the sample solution used for analysis.